




Mahe Medical USATM TRACTA-SCREW

FOR BROKEN & STRIPPED SCREW REMOVAL

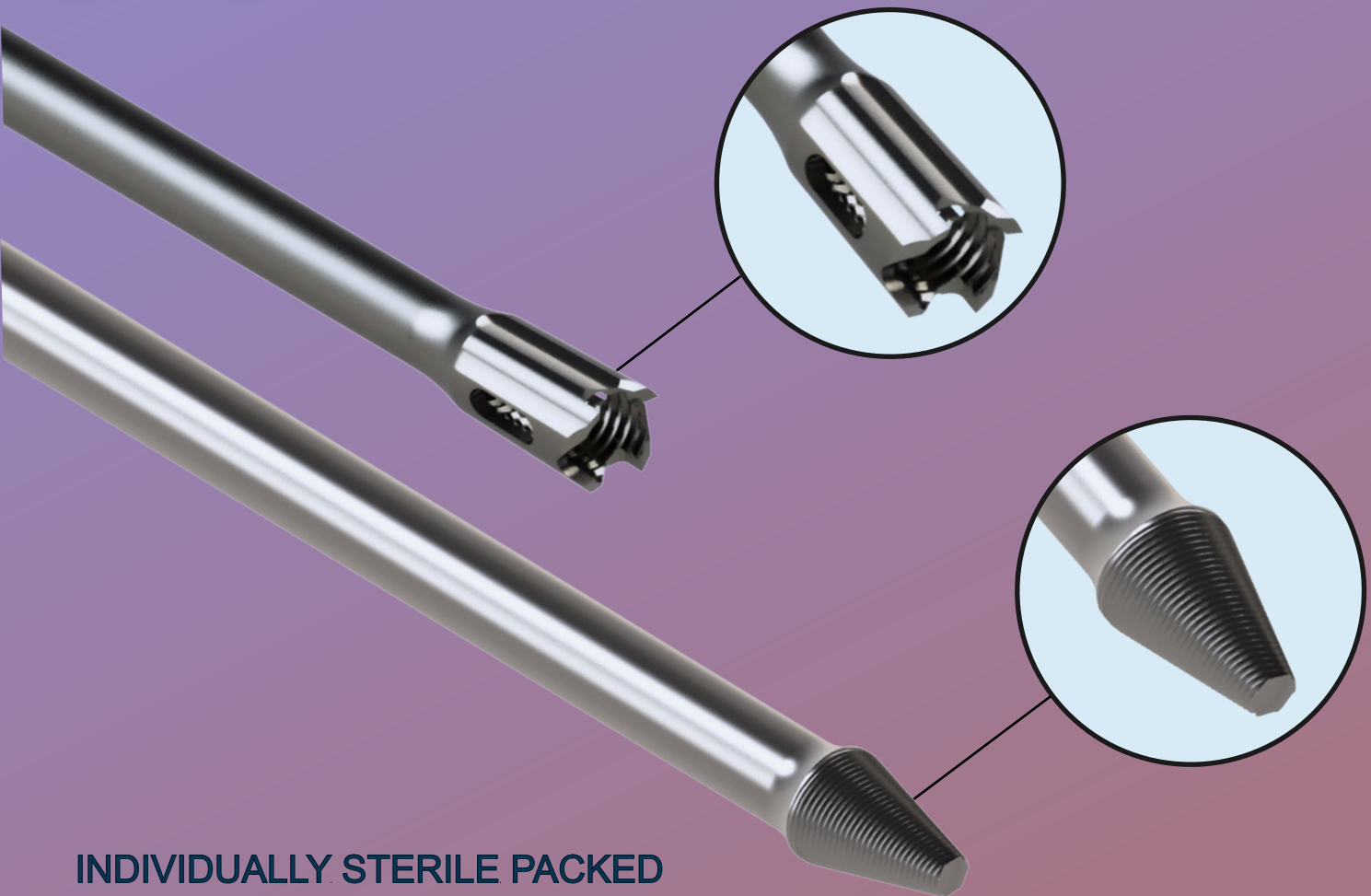

INDIVIDUALLY STERILE PACKED
INNOVATIVE EASY OUT DESIGN

AO QUICK COUPLING



Mahe Medical USATM



Mahe Medical USATM


**EVERY OPERATING ROOM NEEDS
THE SINGLE-USE, STERILE PACKED OPTION
FOR BROKEN & STRIPPED SCREWS
SCREW DIAMETERS 2.0mm - 8.5mm**

FOR BROKEN SCREW REMOVAL

REVERSE ONLY
LEFT TURN
COUNTERCLOCKWISE




Mahe Medical USA[®]
TRACTA-SCREW[®]

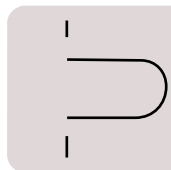

Mahe Medical USA[®]
FOR BROKEN SCREWS
SINGLE USE ONLY


AO QUICK COUPLING

An extractor designed to remove broken screws, with an extended neck length to help reach screws that are deeply buried. Bone clearing windows, reverse cutting threads, and trephine teeth all come together to make the removal quick and to minimize bone loss.

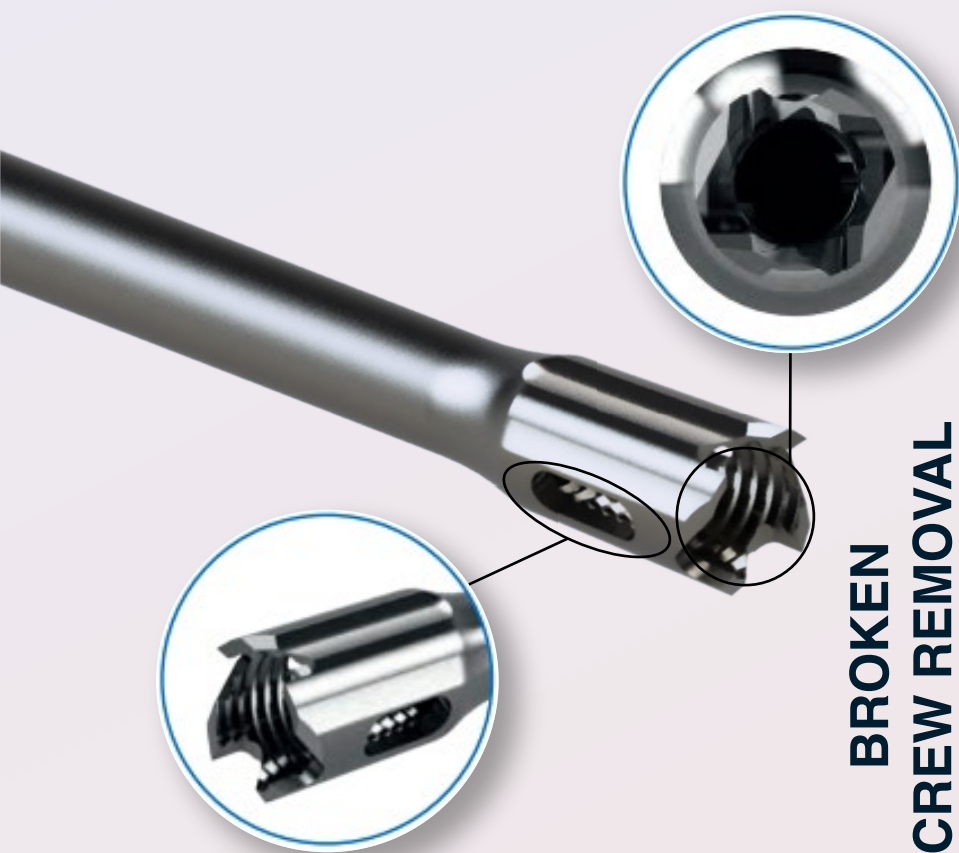

Tapered reverse threads capture the screw quickly while trephine teeth cut through bone.

**BROKEN
SCREW REMOVAL**

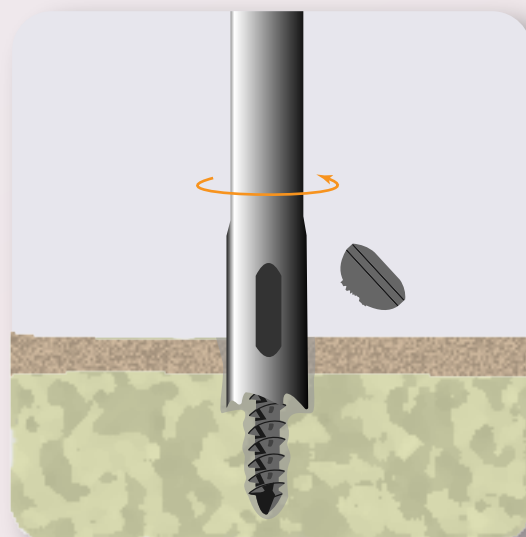


Bone clearing window to allow bone to pass through to prevent blockage.

Single-Use Only: Always use new extractors and drill bits in every procedure. Discard any used single use parts at the conclusion of the case.

FOR STRIPPED SCREW REMOVAL




An extractor designed to remove stripped screws, with a low profile angle to help reach in further. It can even go through spinal tulips to get to the inside of the screw. One-step removal. Simply screw in and turn counterclockwise to remove the screws.

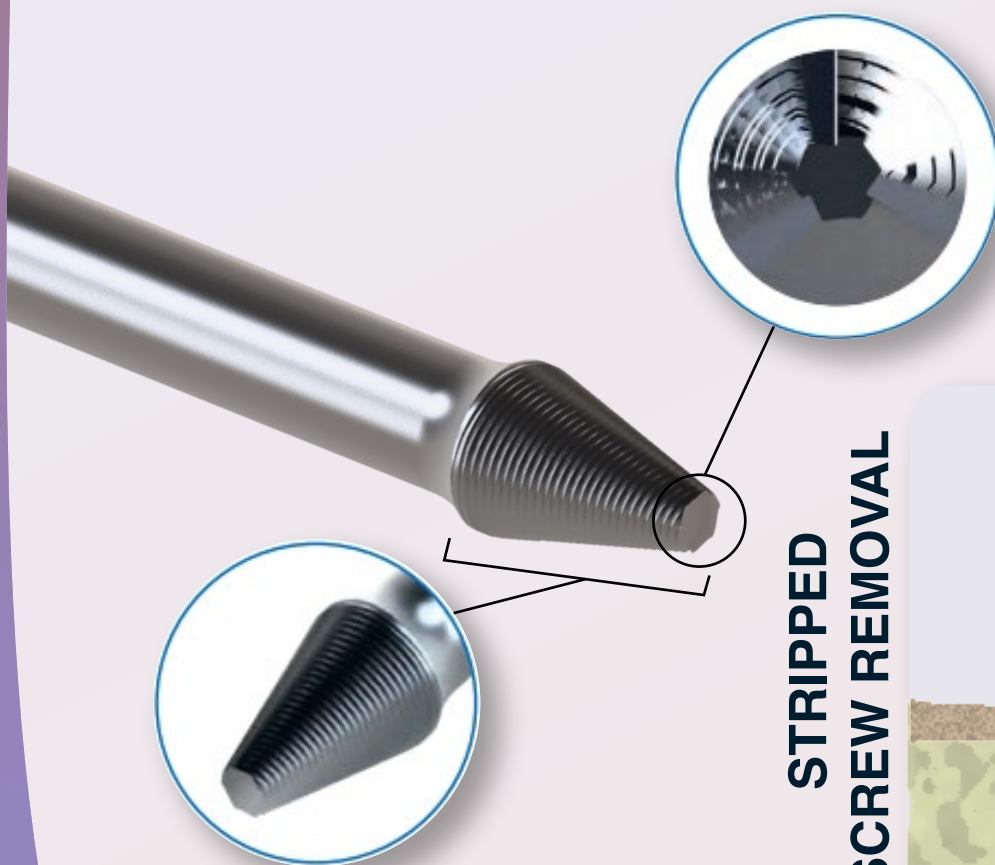

*One-step removal.
Simply screw in and
remove the screw.*

*Reverse threading cuts new
threads into the wall of the
stripped screw well.*

**STRIPPED
SCREW REMOVAL**

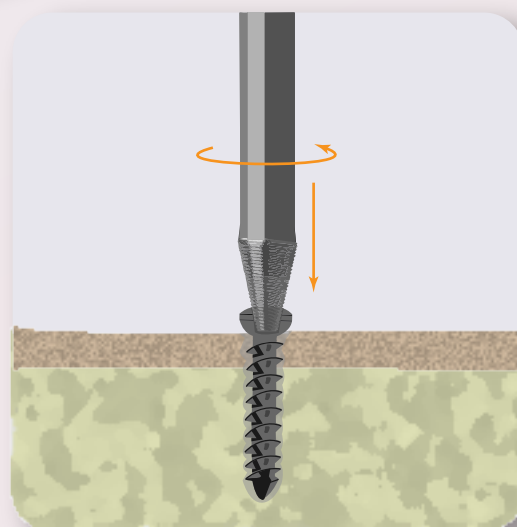








Single-Use Only: Always use new extractors and drill bits in every procedure. Discard any used single use parts at the conclusion of the case.

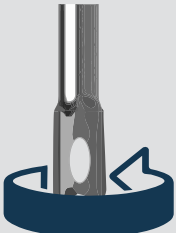








Mahe Medical USATM TRACTA-SCREW

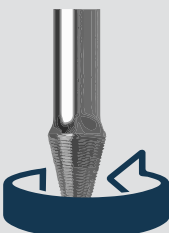FOR BROKEN & STRIPPED SCREW REMOVAL

	FOR BROKEN SCREWS (mm) OUTSIDE THREAD DIAMETER	
	2.0 - 3.0 mm Diameter	
	3.0 - 4.0 mm Diameter	
	4.0 - 5.0 mm Diameter	
	5.0 - 6.0 mm Diameter	
	6.0 - 7.0 mm Diameter	
	7.0 - 8.5 mm Diameter	

Initially select smaller size to begin.

	FOR STRIPPED SCREWS (mm) OUTSIDE THREAD DIAMETER	
	2.0 mm Diameter	
	2.5 - 3.5 mm Diameter	
	3.5 - 4.5 mm Diameter	
	4.5 - 5.5 mm Diameter	
	5.5 - 6.5 mm Diameter	
	6.0 - 7.0 mm Diameter	

Initially select smaller size to begin.

**Xtracta-Screw Sterile Packed
Single-Use Screw Removal Instruments**

Mahe Medical USATM

For Further Information
or to Order, Please Contact:

**3821 Beech St.
Laramie, WY 82070**

**Toll Free: 1 (855) MAHEMED
1 (855) 624-3633**